


# STOP – DO THESE STEPS FIRST

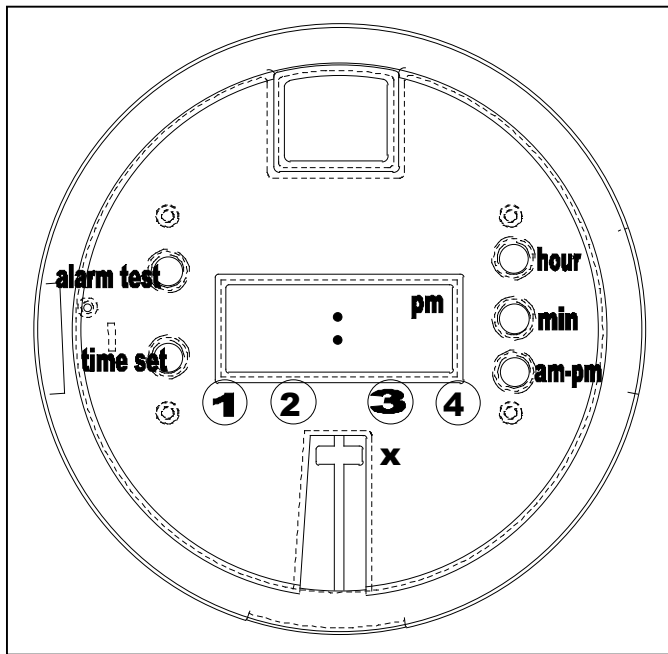

## MedReady 1700PRN Instructions

**1-Install the Battery** The battery compartment is located on the underside of the unit. Open the battery cover by sliding the tab back and lifting after the tab clears the stop. Plug in the battery. Notice the design of the plug allows it to be inserted in only one direction. Replace the battery cover. If the battery has been installed correctly the colon on the display face should now be flashing. \*\*\* **The battery must be connected for the MedReady to operate** \*\*\*

**2 – AC Adapter** Plug the AC adapter into the underside of the unit where it says “12 Volt DC Jack Here” and connect to an active electrical outlet, which is not controlled by a wall switch. The entire LED display should now be visible.

**3 – Remove the Blue Lid** Insert the round key into the lock, depress slightly and rotate the key

clockwise to the open / horizontal position. With the MedReady in front of you, place one hand at the 9 o'clock position and one hand at the 3 o'clock position. Depress slightly with both hands and **rotate clockwise** about one inch until the blue top stops rotating and then lift the blue cover off. Leave the key in the blue top, and in the open / horizontal position.

**4 – Setting the PRN Time** While continuing to hold down the ‘1’ button, adjust the PRN time by pressing the HOUR or MIN until the desired PRN time is displayed. Note that only one PRN time can be set, therefore the 2, 3 and 4 buttons have been disabled.

**5 – Loading Medication** The medication tray can be removed by simply lifting it out of the dispenser. There are 29 medication positions of which 28 are usable.

Starting with any empty slot, load each slot consecutively in a counterclockwise direction until all 28 usable positions are filled, or until you have reached your desired fill quantity. Do not skip any slots while loading the medication. It is alright to leave empty slots at the end. Example: If loading only 24 doses, then 24 doses should be loaded in a counterclockwise direction leaving the remainder of the slots empty.

It important to position the medication tray with the empty position at the 6:00 position and centered directly under the white T-Bar.

**IF LOADING LESS THAN 28 Doses:** Since the medication tray rotates clockwise as it dispenses, if all 28 doses are not loaded, the medication tray must be positioned with the last empty slot directly under the white T-Bar, and medication starting directly to the right (counterclockwise).

## 6 – Replacing the Blue Lid.

**Before putting the blue top on, look at the Red Cover and make sure the white T-BAR is pushed up towards the arrow “→”.** Insure the round key is in the lock and in the open / horizontal position. **Close the Medication door on the Blue Lid.** Place the blue lid on the MedReady slightly rotated clockwise about an inch (in the same position in which it was removed). With the MedReady in front of you, place one hand at the 9 O'clock position and one hand at the 3 O'clock position. Depress **slightly** with both hands and rotate **counter clockwise** about one inch until the blue top tabs slide under the base tabs. Rotate the key **counter clockwise** to the closed / vertical position and remove.

## Other Settings and Functions:

- **Away Mode** You may turn the dispenser **OFF** by pressing the *hour* button and holding for approximately 4 seconds. The display will read '**OFF**' and the medication will stop dispensing. To return to normal operation, press the *hour* until the PRN time displays.

- **Alarm Test Press** and release this button on the upper left to test the dispensing rotation. Slide the **white T-BAR** up reset the PRN time.

## Trouble Shooting

No PRN time is being displayed, but the colon is blinking	The Device is not getting A/C power. Try another A/C wall outlet that you know is live. Insure both ends are connected.
The PRN Time does not reset when opening the medication door.	Instruction # 6 was not followed when replacing the blue lid. The medication door pushes up the white T-BAR (directly below the medication door), which triggers a micro switch and resets the PRN time. If the medication door is on the wrong side of the white T-BAR, the medication door can't engage the micro switch to reset the PRN time.
The medication tray does not rotate when the PRN time reaches 0:00	The battery is not plugged in. See step # 1. To test the battery connection, disconnect the A/C power, and <b>ONLY</b> the colon should be blinking.
The new medication slot is empty	The medication was not loaded properly, Refer to step #5.
Display shows 'OFF'	Depress the Hour button by itself for several seconds to return to show the PRN time.
Blue Lid won't go back on	See # 6 above
The Battery and A/C are connected but the device won't rotate to the next medication when the PRN time reaches 0:00	Most likely, the medication tray was previously forced to the next dose disabling the device. Please call our customer support at 310.328.7557 (9-4:30 M-F Pacific Time) to rectify this situation.

## Please read and understand the following:

Note that many PRN regimens involve the dispensing of narcotics. The MedReady 1700PRN is designed as a medication aid only. MedReady is NOT a medication Safe or Vault. With determination, a patient will be able to break into the device and access more than intended. If there is any doubt that a patient will have an uncontrollable urge to access more medication than prescribed, please do NOT use this device. Under normal use, only one dose is available at any one time.